

# Capillary Line Fixture™

## Capillary Line Fixture™

Model ECT/CLF/A

### Main Features:

- Center-of-rotation (COR) error measurement for systems that have axially dependent COR correction
- Evaluation of changes of radius-of-rotation on spatial resolution
- Planar and SPECT spatial resolution measurements in air
- Quantitative evaluation of reconstruction filters
- Research


### Specifications:

Overall length of line fixture 55.9 cm

Usable length 43.2 cm

Outer support tubing: Polycarbonate

Inner diameter of tubing: 1.2 mm

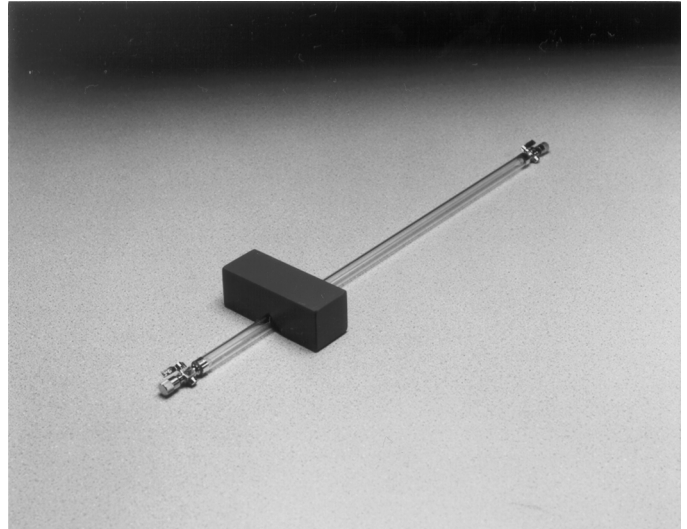

*Capillary Line Fixture™*

### Also Available:

## Capillary Line Fixture, Special Length™

Model ECT/CLF-X/A

Same as above, but length cut to order